


**This is TERUMO Territory.™**

## TERUMO *Glidewire*®

**What makes Terumo Glidewire perform like the Glidewire you know and trust?**

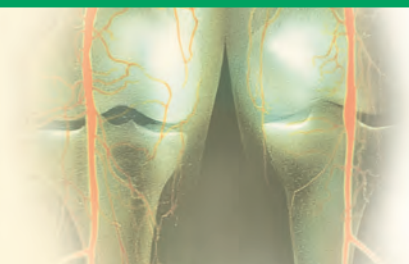

### Terumo Glide Technology™ hydrophilic coating

- Unmatched lubricity consistently gets you where you need to go

### Core-to-tip design

- Provides optimal torque control for superior navigation


### Polyurethane jacket with tungsten

- Provides enhanced radiopacity of entire Terumo Glidewire

### Terumo nitinol alloy core

- Resists kinking for easier, faster navigation through even the most tortuous vessels
- Restores shape to manufactured state for multiple use in a single procedure

**Consistent, reliable Terumo Glidewire performance that gets you there**

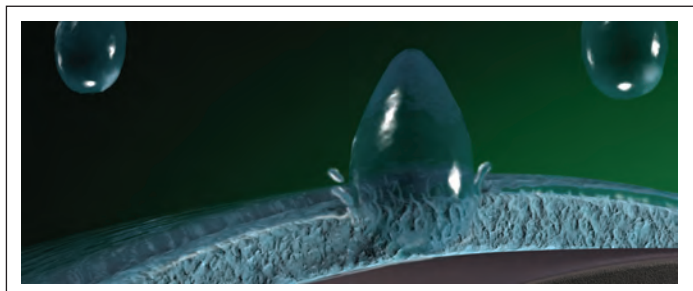

### Terumo Glide Technology hydrophilic coating

- Unique hydrophilic coating has been imitated but never equaled
- Consistent performance gets you there faster

## Peripheral Guidewires (0.032" - 0.038")

# Terumo *Glidewire*<sup>®</sup>

### Hydrophilic Coated Guidewires

Packed 5 per box.

	Product Code	Wire Type	Diameter (in)	Total Length (cm)	Flexible Tip Length (cm)	Tip Shape
<b>Standard Glidewire</b>	GR3202	Standard	0.032	150	3	Angle
	GR3201	Standard	0.032	150	3	Straight
	GR3513	Standard	0.035	80	3	Angle
	GR3505	Standard	0.035	120	3	Angle
	GR3506	Standard	0.035	150	3	Angle
	GR3508	Standard	0.035	180	3	Angle
	GR3509	Standard	0.035	260	3	Angle
	GR3501	Standard	0.035	150	3	Straight
	GR3503	Standard	0.035	180	3	Straight
	GR3504	Standard	0.035	260	3	Straight
	GR3811	Standard	0.038	80	3	Angle
	GR3805	Standard	0.038	120	3	Angle
	GR3806	Standard	0.038	150	3	Angle
	GR3808	Standard	0.038	180	3	Angle
	GR3809	Standard	0.038	260	3	Angle
	GR3801	Standard	0.038	150	3	Straight
	GR3803	Standard	0.038	180	3	Straight
	GR3804	Standard	0.038	260	3	Straight
<b>Shapeable Tip Glidewire</b>	GE3501	Shapeable	0.035	150	3	Straight/Shapeable
	GE3502	Shapeable	0.035	180	3	Straight/Shapeable
	GE3801	Shapeable	0.038	150	3	Straight/Shapeable
<b>Long Taper Glidewire</b>	GR3507	Long Taper	0.035	150	8	Angle
	GR3502	Long Taper	0.035	150	8	Straight
	GR3807	Long Taper	0.038	150	8	Angle
	GR3802	Long Taper	0.038	150	8	Straight
<b>1 cm Taper Glidewire</b>	GR3515	Standard	0.035	150	1	Angle
	GR3516	Standard	0.035	180	1	Angle
	GR3517	Standard	0.035	260	1	Angle
	GR3518	Standard	0.035	150	1	Straight
	GR3519	Standard	0.035	180	1	Straight
	GR3520	Standard	0.035	260	1	Straight
	GS3511	Stiff	0.035	150	1	Angle

Peripheral Guidewires (0.032" – 0.038")

# Terumo *Glidewire*®

**Hydrophilic Coated Guidewires**

Packed 5 per box.

	Product Code	Wire Type	Diameter (in)	Total Length (cm)	Flexible Tip Length (cm)	Tip Shape
<b>Stiff Shaft Glidewire</b>	GS3505	Stiff	0.035	80	3	Angle
	GS3506	Stiff	0.035	150	3	Angle
	GS3508	Stiff	0.035	180	3	Angle
	GS3509	Stiff	0.035	260	3	Angle
	GS3501	Stiff	0.035	150	3	Straight
	GS3503	Stiff	0.035	180	3	Straight
	GS3504	Stiff	0.035	260	3	Straight
	GS3804	Stiff	0.038	80	3	Angle
	GS3805	Stiff	0.038	150	3	Angle
	GS3807	Stiff	0.038	180	3	Angle
	GS3801	Stiff	0.038	150	3	Straight
	GS3803	Stiff	0.038	180	3	Straight
	GS3511	Stiff	0.035	150	1	Angle
<b>Stiff Shaft Long Taper Glidewire</b>	GS3507	Stiff	0.035	150	5	Angle
	GS3510	Stiff	0.035	260	5	Angle
	GS3502	Stiff	0.035	150	5	Straight
	GS3806	Stiff	0.038	150	5	Angle
	GS3802	Stiff	0.038	150	5	Straight
<b>J-Tip Glidewire</b>	GR3510	Standard	0.035	150	3	J-Tip
	GR3810	Standard	0.038	150	3	J-Tip
<b>Bolia Curve Glidewire</b>	GR3521	Standard	0.035	180	5	Bolia
	GR3522	Standard	0.035	180	8	Bolia

## Small Vessel Guidewires (0.018" – 0.025")

# Terumo *Glidewire*®

### Hydrophilic Coated Guidewires

Packed 5 per box.

	Product Code	Wire Type	Diameter (in)	Total Length (cm)	Flexible Tip Length (cm)	Tip Shape
<b>Standard Glidewire</b>	GR1804	Standard	0.018	80	3	Angle
	GR1805	Standard	0.018	150	3	Angle
	GR1806	Standard	0.018	180	3	Angle
	GR1802	Standard	0.018	150	3	Straight
	GR1803	Standard	0.018	180	3	Straight
	GR1807	Standard	0.018	260	3	Angle
	GR1808	Standard	0.018	260	3	Straight
	GR2504	Standard	0.025	150	3	Angle
	GR2505	Standard	0.025	180	3	Angle
	GR2506	Standard	0.025	260	3	Angle
	GR2501	Standard	0.025	150	3	Straight
	GR2502	Standard	0.025	180	3	Straight
	GR2503	Standard	0.025	260	3	Straight
<b>Shapeable Tip Glidewire</b>	GE1801	Shapeable	0.018	150	3	Straight/Shapeable

## Terumo – A Global Leader in Medical Devices

### Terumo Introducers

Pinnacle® Introducer Sheaths  
 Pinnacle® Destination® Guiding Sheaths  
 Pinnacle® R/O II Radiopaque Marker Introducer Sheaths  
 Pinnacle® R/O II HiFlo Introducer Sheaths  
 Pinnacle TIF Tip™ Introducer Sheaths  
 Glidesheath™ Introducer Sheaths  
 GlideAccess™ Micro Accessing Systems

### Terumo Guidewires

Glidewire® Hydrophilic Coated Guidewires  
 Glidewire® GT Hydrophilic Coated Guidewires  
 Glidewire® Gold Hydrophilic Coated Guidewires  
 Runthrough™ NS Coronary Guidewire

### Terumo Angiographic Catheters

Glidecath® Hydrophilic Coated Catheters  
 Glidecath® XP Hydrophilic Coated Catheters

### Terumo Microcatheters

Progreat™ Hydrophilic Coated Microcatheters  
 Progreat™ Coaxial Microcatheter Systems

### Terumo Radial Compression Devices

TR Band™ Radial Compression Devices

### Terumo Interventional Products

AZUR™ Peripheral HydroCoil® Embolization Systems  
 Bead Block™ Compressible Microspheres

For more information, call 800.862.4143 to speak to an Inside Sales Customer Care Representative, or visit [www.terumo.is.com](http://www.terumo.is.com)  
 Terumo Interventional Systems · 2101 Cottontail Lane · Somerset, NJ 08873 · Fax: 800.411.5870

©2008 Terumo Medical Corporation. All rights reserved. HydroCoil is a registered trademark of MicroVenton Inc. Bead Block is a trademark of Biocompatibles UK Ltd.  
 All other brand names are trademarks or registered trademarks of Terumo. Lit. # TIS003 05/08