

*A master in the art of
multipurpose microcatheters.*

MERIT
Maestro® MICROCATHETER



MERIT MAESTRO[®]

MICROCATHETER

*A Multipurpose Microcatheter
designed for use in peripheral
and coronary vasculature.*

FEATURE & BENEFITS

FLEXIBLE DISTAL REGION 20 CM

designed for atraumatic vessel entry

HYDROPHILIC COATING DISTAL 80 CM

for smooth traction through tortuous anatomy

RADIOPAQUE MARKER

for easy identification of distal tip of catheter.

Located 1.3 mm proximal to the distal tip

NEW SWAN NECK DESIGN

high performance of preshaped distal end

ideal for access of hepatic artery




(2) 3 ml Medallion[®] syringes included for ease of mixing embolics/spheres



NEW

**MERIT MEDICAL'S
SWAN NECK DESIGN**

*Manufactured
pre-shaped
Swan Neck design.*

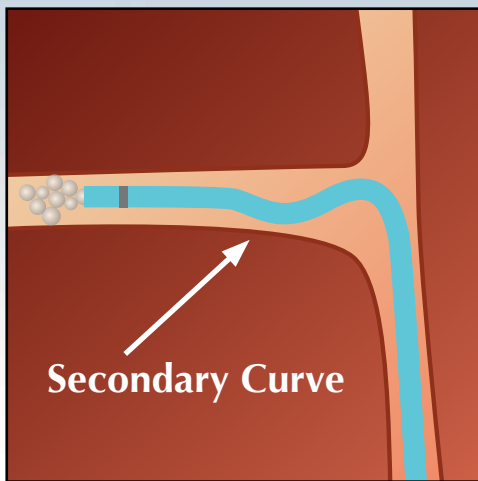


*Once the Swan Neck
catheter tip is in
position in the vessel,
the secondary curve
helps seat the catheter
in the vessel reducing
the recoiling effect of
the embolic agent as
it is introduced.*

Merit Maestro[®]

MICROCATHETER

ORDERING INFORMATION

Catalog Number	Proximal	Distal	Length cm	Tip	Inner Diameter	Guide Wire	Guide Catheter Indicated
2 . 8 F T A P E R E D T O 2 . 4 F							
28MC24110ST	2.8F	2.4F	110	Straight	.020" (0.52 mm)	.018" (0.46 mm)	.040" (1.02 mm)
28MC24130ST	2.8F	2.4F	130	Straight	.020" (0.52 mm)	.018" (0.46 mm)	.040" (1.02 mm)
28MC24150ST	2.8F	2.4F	150	Straight	.020" (0.52 mm)	.018" (0.46 mm)	.040" (1.02 mm)
28MC2411045	2.8F	2.4F	110	45 degree	.020" (0.52 mm)	.018" (0.46 mm)	.040" (1.02 mm)
28MC2413045	2.8F	2.4F	130	45 degree	.020" (0.52 mm)	.018" (0.46 mm)	.040" (1.02 mm)
28MC2415045	2.8F	2.4F	150	45 degree	.020" (0.52 mm)	.018" (0.46 mm)	.040" (1.02 mm)
28MC24110SN	2.8F	2.4F	110	Swan Neck	.020" (0.52 mm)	.018" (0.46 mm)	.040" (1.02 mm)
28MC24130SN	2.8F	2.4F	130	Swan Neck	.020" (0.52 mm)	.018" (0.46 mm)	.040" (1.02 mm)
28MC24150SN	2.8F	2.4F	150	Swan Neck	.020" (0.52 mm)	.018" (0.46 mm)	.040" (1.02 mm)
2 . 8 F N O T A P E R							
28MC28110ST	2.8F	2.8F	110	Straight	.024" (0.62 mm)	.021" (0.53 mm)	.040" (1.02 mm)
28MC28130ST	2.8F	2.8F	130	Straight	.024" (0.62 mm)	.021" (0.53 mm)	.040" (1.02 mm)
28MC28150ST	2.8F	2.8F	150	Straight	.024" (0.62 mm)	.021" (0.53 mm)	.040" (1.02 mm)
28MC2811045	2.8F	2.8F	110	45 degree	.024" (0.62 mm)	.021" (0.53 mm)	.040" (1.02 mm)
28MC2813045	2.8F	2.8F	130	45 degree	.024" (0.62 mm)	.021" (0.53 mm)	.040" (1.02 mm)
28MC2815045	2.8F	2.8F	150	45 degree	.024" (0.62 mm)	.021" (0.53 mm)	.040" (1.02 mm)
28MC28110SN	2.8F	2.8F	110	Swan Neck	.024" (0.62 mm)	.021" (0.53 mm)	.040" (1.02 mm)
28MC28130SN	2.8F	2.8F	130	Swan Neck	.024" (0.62 mm)	.021" (0.53 mm)	.040" (1.02 mm)
28MC28150SN	2.8F	2.8F	150	Swan Neck	.024" (0.62 mm)	.021" (0.53 mm)	.040" (1.02 mm)
2 . 9 F N O T A P E R							
29MC29110ST	2.9F	2.9F	110	Straight	.027" (0.68 mm)	.021" (0.53 mm)	.042" (1.07 mm)
29MC29130ST	2.9F	2.9F	130	Straight	.027" (0.68 mm)	.021" (0.53 mm)	.042" (1.07 mm)
29MC29150ST	2.9F	2.9F	150	Straight	.027" (0.68 mm)	.021" (0.53 mm)	.042" (1.07 mm)
29MC2911045	2.9F	2.9F	110	45 degree	.027" (0.68 mm)	.021" (0.53 mm)	.042" (1.07 mm)
29MC2913045	2.9F	2.9F	130	45 degree	.027" (0.68 mm)	.021" (0.53 mm)	.042" (1.07 mm)
29MC2915045	2.9F	2.9F	150	45 degree	.027" (0.68 mm)	.021" (0.53 mm)	.042" (1.07 mm)
29MC29110SN	2.9F	2.9F	110	Swan Neck	.027" (0.68 mm)	.021" (0.53 mm)	.042" (1.07 mm)
29MC29130SN	2.9F	2.9F	130	Swan Neck	.027" (0.68 mm)	.021" (0.53 mm)	.042" (1.07 mm)
29MC29150SN	2.9F	2.9F	150	Swan Neck	.027" (0.68 mm)	.021" (0.53 mm)	.042" (1.07 mm)

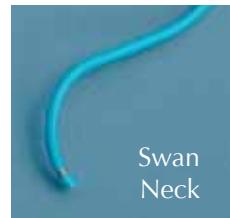

All Merit Maestro Microcatheters are rated to 800psi infusion dynamic pressure.

	EMBOLICS		Coils
	Particles	Spherical	
2.8F / 2.4F	≤ 700 um Emboli	≤ 700 um Microspheres	0.46mm / 0.018" 0.46mm / 0.018"
2.8F / 2.8F	≤ 700 um Emboli	≤ 700 um Microspheres	
2.9F / 2.9F	≤ 900 um Emboli	≤ 900 um Microspheres	

