

EN Snare®

Endovascular Snare System

Capture the Possibilities™


EN Snare®

The EN Snare Endovascular System is intended for use in the cardiovascular system or hollow viscous to retrieve and manipulate foreign objects.

Manipulation procedures include indwelling venous catheter repositioning, indwelling venous catheter fibrin sheath stripping and central venous access venipuncture procedure assistance.


FEATURES

- Available in seven sizes with overlapping diameters to accommodate a broad range of vessels throughout the body.
- The EN Snare is designed to rotate and expand in confined spaces to retrieve objects.
- The three preformed loops of the EN Snare are interlaced and self-sizing for excellent vessel coverage.
- Super-elastic Nitinol wire provides flexibility, kink-resistance, and 1:1 torque control.
- Interlaced Nitinol loops provide structural integrity.
- Platinum strands are incorporated in the snare loops for excellent fluoroscopic visualization.

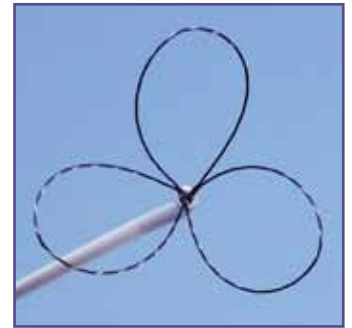


Snare Catheter:


FEP, kink-resistant guiding catheter includes radiopaque band that is easily visible during fluoroscopy. The 6 French and 7 French snare catheters are designed with a unique 15 degree angled tip, allowing the device to be steered during placement.

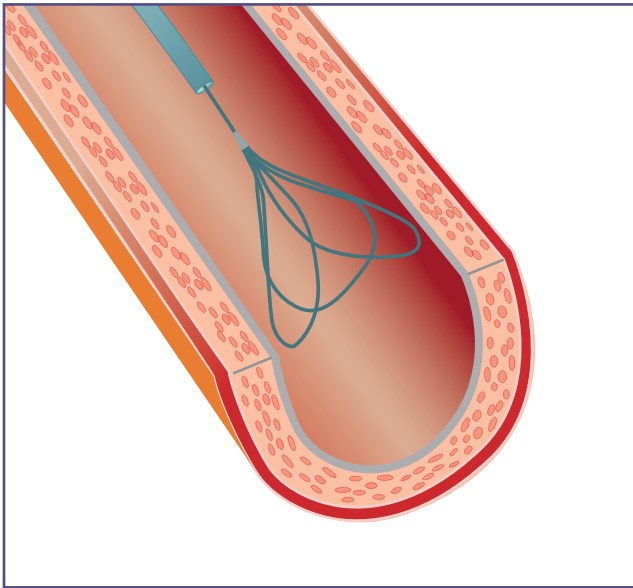

The EN Snare is designed to expand to vessel wall size for retrieval of hard to get side-walled objects.

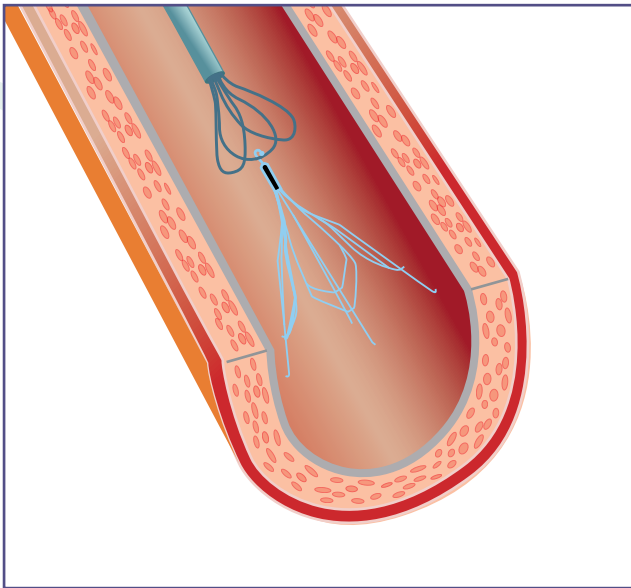

The EN Snare is designed with three loops to increase the probability of foreign body capture.

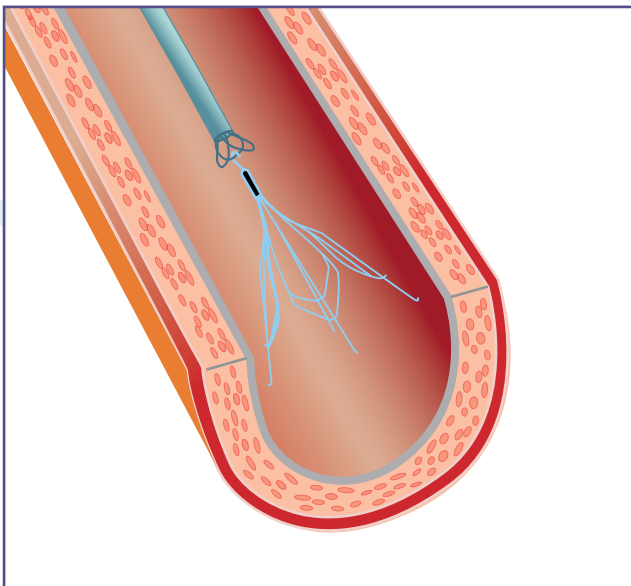

Interlaced Nitinol loops provide structural integrity when capturing objects.



ORDERING INFORMATION

Catalog Number	Working Diameter	Snare Length	French Size	Snare Catheter Length
EN1003004 Mini	2-4mm	175cm	3.2Fr	150cm
EN1003008 Mini	4-8mm	175cm	3.2Fr	150cm
EN2006010 Standard	6-10mm	120cm	6Fr	100cm
EN2006015 Standard	9-15mm	120cm	6Fr	100cm
EN2006020 Standard	12-20mm	120cm	6Fr	100cm
EN2007030 Standard	18-30mm	120cm	7Fr	100cm
EN2007045 Standard	27-45mm	120cm	7Fr	100cm

The **EN Snare System** is sold individually and includes a snare, torque device, introducer and snare catheter.

Snare and catheter are packaged in separate pouches.



Merit Medical Systems, Inc. • 1600 West Merit Parkway • South Jordan, Utah 84095 • 1-801-253-1600 • 1-800-35-MERIT

www.merit.com

Merit Medical Europe • Horsterweg 24 • 6199 AC MAASTRICHT- AIRPORT The Netherlands • +31 43 35 88 222

Merit Medical Ireland Ltd. • Parkmore Business Park West Galway, Ireland • +353 (0) 91 703 733

Free phone from specific country:

Belgium 0800 72906 • France 0800 91 60 30 • Germany 0800 182 0871 • Luxembourg 0800 22522 • Netherlands 0800 022 8184 • UK 0800 973 115

401896001/A ID 113009