



FOUNTAIN® INFUSION SYSTEMS

Better by Design

GRADIENT HOLE SIZES

Exclusive gradient hole sizes allow uniform infusion of therapeutic agent along the entire length of the infusion segment. Holes are laser drilled to provide a superior infusion pattern for continuous slow-drip or pulsed infusion techniques.

LONGER INFUSION SEGMENTS

Both the 4F and 5F products have infusion segments from 5cm - 50cm.

CATHETER

Compatible with 0.035" (0.89mm) guide wire. Fabricated to provide the proper stiffness for good trackability and ease of placement.

OCCUDING WIRE



Effectively blocks off the end hole of the catheter so the therapeutic agent is directed out the side holes into the thrombus.

HEMOSTASIS VALVE

Large bore valve allows easy wire placement while minimizing patient blood loss. Airless rotator minimizes bubbles.

HIGH PRESSURE

Check Relief Valve
Designed so injection occurs only from the infusion syringe. Unique "in-line" design provides better handling during pulsed infusion.

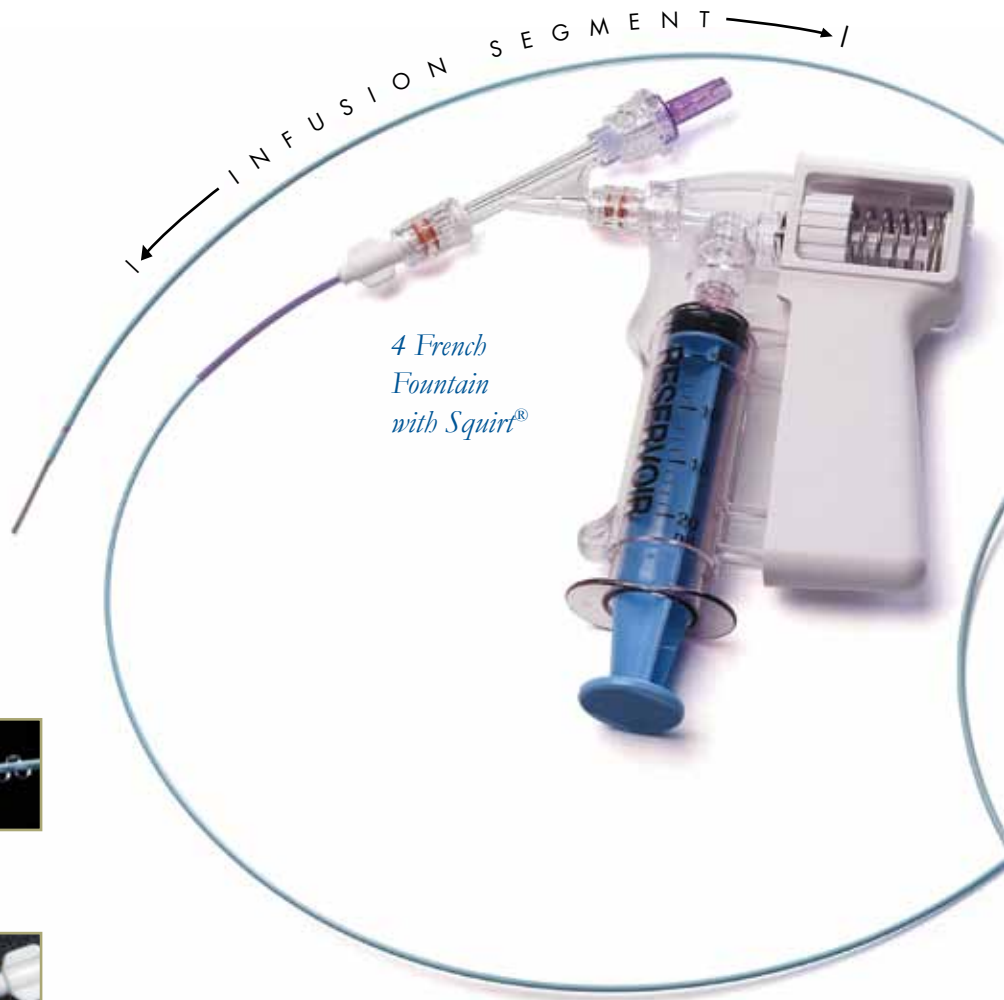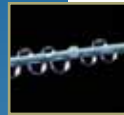

4F PROPRIETARY SNAP CAP

The cap at the proximal end of the occluding wire snaps into the hemostasis valve ensuring correct wire placement and adding security for overnight or long-term infusions.

FOUNTAIN[®]

INFUSION SYSTEMS

*4 French
Fountain
with Squirt[®]*

THE 4 FRENCH AND 5 FRENCH FOUNTAIN INFUSION SYSTEMS ARE AVAILABLE WITH THE SQUIRT OR WITH AN INFUSION SYSTEM CONTAINING A HIGH PRESSURE CHECK-RELIEF VALVE AND INFUSION SYRINGE. BOTH SYSTEMS CONTAIN A MERIT ACCESSPLUS[™] HEMOSTASIS VALVE AND A 20ML RESERVOIR SYRINGE.




5F PROPRIETARY WIRE PROTECTOR CAP

The wire protector cap covers the proximal end of the occluding wire and snaps into the hemostasis valve.

The side holes on the Fountain Infusion Catheter rotate around the catheter to create a more uniform dispersion of therapeutic agent.

FOUNTAIN

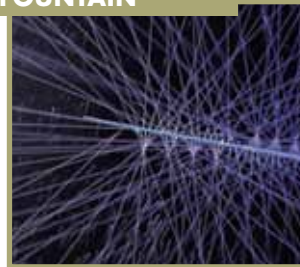

COMPETITION



Other competitive catheters have side holes or slits lined up in 4 parallel rows



SQUIRT[®] FLUID DELIVERY SYSTEM

Consistent, forceful, pulsed injections for optimal thrombolysis procedures.

Ordering Information

FOUNTAIN INFUSION SYSTEM WITH SQUIRT® FLUID DELIVERY SYSTEM

Catalog Number	French Size	Catheter Length cm	Infusion Segment cm
FIS4-45-10SQ	4F	45	10
FIS4-45-20SQ	4F	45	20
FIS4-90-5SQ	4F	90	5
FIS4-90-10SQ	4F	90	10
FIS4-90-20SQ	4F	90	20
FIS4-90-30SQ	4F	90	30
FIS4-90-40SQ	4F	90	40
FIS4-90-50SQ	4F	90	50
FIS4-135-5SQ	4F	135	5
FIS4-135-10SQ	4F	135	10
FIS4-135-20SQ	4F	135	20
FIS4-135-30SQ	4F	135	30
FIS4-135-40SQ	4F	135	40
FIS4-135-50SQ	4F	135	50

FIS5-45-10SQ	5F	45	10
FIS5-45-20SQ	5F	45	20
FIS5-90-5SQ	5F	90	5
FIS5-90-10SQ	5F	90	10
FIS5-90-20SQ	5F	90	20
FIS5-90-30SQ	5F	90	30
FIS5-90-40SQ	5F	90	40
FIS5-90-50SQ	5F	90	50
FIS5-135-5SQ	5F	135	5
FIS5-135-10SQ	5F	135	10
FIS5-135-20SQ	5F	135	20
FIS5-135-30SQ	5F	135	30
FIS5-135-40SQ	5F	135	40
FIS5-135-50SQ	5F	135	50

Kit includes Fountain catheter, occluding wire, AccessPLUS hemostasis valve, Squirt Fluid Delivery Syringe and 20ml reservoir syringe. (wire protector cap with 5F)

FOUNTAIN INFUSION SYSTEM

Catalog Number	French Size	Catheter Length cm	Infusion Segment cm
FIS4-45-10	4F	45	10
FIS4-45-20	4F	45	20
FIS4-90-5	4F	90	5
FIS4-90-10	4F	90	10
FIS4-90-20	4F	90	20
FIS4-90-30	4F	90	30
FIS4-90-40	4F	90	40
FIS4-90-50	4F	90	50
FIS4-135-5	4F	135	5
FIS4-135-10	4F	135	10
FIS4-135-20	4F	135	20
FIS4-135-30	4F	135	30
FIS4-135-40	4F	135	40
FIS4-135-50	4F	135	50

FIS5-45-10	5F	45	10
FIS5-45-20	5F	45	20
FIS5-90-5	5F	90	5
FIS5-90-10	5F	90	10
FIS5-90-20	5F	90	20
FIS5-90-30	5F	90	30
FIS5-90-40	5F	90	40
FIS5-90-50	5F	90	50
FIS5-135-5	5F	135	5
FIS5-135-10	5F	135	10
FIS5-135-20	5F	135	20
FIS5-135-30	5F	135	30
FIS5-135-40	5F	135	40
FIS5-135-50	5F	135	50

Kit includes Fountain catheter, occluding wire, AccessPLUS hemostasis valve, high pressure check-relief valve, 1ml infusion syringe and 20ml reservoir syringe. (wire protector cap with 5F)

SQUIRT FLUID DELIVERY SYRINGE

Catalog Number	FDS-100
----------------	---------

CATHETER, OCCLUDING WIRE AND AccessPLUS™ HEMOSTASIS VALVE

Catalog Number	Size	Catheter French cm	Infusion Length cm	Segment
IS4-45-10		4F	45	10
IS4-45-20		4F	45	20
IS4-90-5 4F		90	5	
IS4-90-10		4F	90	10
IS4-90-20		4F	90	20
IS4-90-30		4F	90	30
IS4-90-40		4F	90	40
IS4-90-50		4F	90	50
IS4-135-5		4F	135	5
IS4-135-10		4F	135	10
IS4-135-20		4F	135	20
IS4-135-30		4F	135	30
IS4-135-40		4F	135	40
IS4-135-50		4F	135	50

IS5-45-10		5F	45	10
IS5-45-20		5F	45	20
IS5-90-5 5F		90	5	
IS5-90-10		5F	90	10
IS5-90-20		5F	90	20
IS5-90-30		5F	90	30
IS5-90-40		5F	90	40
IS5-90-50		5F	90	50
IS5-135-5		5F	135	5
IS5-135-10		5F	135	10
IS5-135-20		5F	135	20
IS5-135-30		5F	135	30
IS5-135-40		5F	135	40
IS5-135-50		5F	135	50

Includes catheter, occluding wire and AccessPLUS hemostasis valve. (wire protector cap with 5F)

OCCLUDING WIRE

Catalog Number	Catheter French Size	Catheter Length cm
FW4-45	4F	45
FW4-90	4F	90
FW4-135	4F	135
FW5-45 5F	45	
FW5-90 5F	90	
FW5-135	5F	135



Merit Medical Systems, Inc. • 1600 West Merit Parkway, South Jordan, Utah 84095 • 1-801-253-1600 • 1-800-35-MERIT
 Merit Medical Europe • Horsterweg 24, 6199 AC MAASTRICHT- AIRPORT, The Netherlands • +31 43 35 88 222
 Merit Medical Ireland Ltd. • Parkmore Business Park West, Galway, Ireland • +353 (0) 91 703 733
 Free phone for specific country: Austria 0800 295 374 • Belgium 0800 72 906 (Dutch) 0800 73 172 (French) • Denmark 80 88 00 24 • France 0800 91 60 30
 Germany 0800 182 0871 • Luxembourg 0800 22 522 • Netherlands 0800 022 81 84 • Sweden 020 792 445 • UK 0800 973 115

www.merit.com